

FemTouch™

on AcuPulse™

360°
of Vaginal Health



FemTouch™:

The 360° solution

for vaginal health-related conditions

A mutual collaboration with leading experts resulted in the development of FemTouch – a vaginal laser treatment driven by practitioners and women’s needs

How will FemTouch benefit your patients?

CO₂ laser treatments have been shown to promote restoration of the physiology of the vaginal canal to a premenopausal-like stage and to the reduction of the vulvo-vaginal ageing symptoms that accompany vaginal atrophy and incontinence.



“My menopause patients are so thankful for helping to restore their sex lives and relieve their vaginal dryness. One patient said she felt 19 again!”

Dr. Scott Eder, GYN, NJ, USA



EFFICIENT

Only 2-4 treatments needed for optimal outcomes



HASSLE-FREE

No downtime or post-procedure care required



FAST

Comfortable IPT and OPT procedure lasting just a few minutes



VERSATILE


Combined with AcuPulse for various aesthetic and gynecological indications

How does FemTouch work?

FemTouch addresses vaginal health related conditions by promoting the remodeling of the vaginal mucosa. The Fractional CO₂ laser is applied in 360° along the vaginal canal which results in a very gentle micro-ablation of the vaginal lining.

The energy level of 7.5mJ - 12.5mJ utilized working with FemTouch ensures that the desirable clinical effect is achieved while the fibromuscular layer remains intact.

Histological effect on vaginal mucosa* - impact depth 200-500 µm

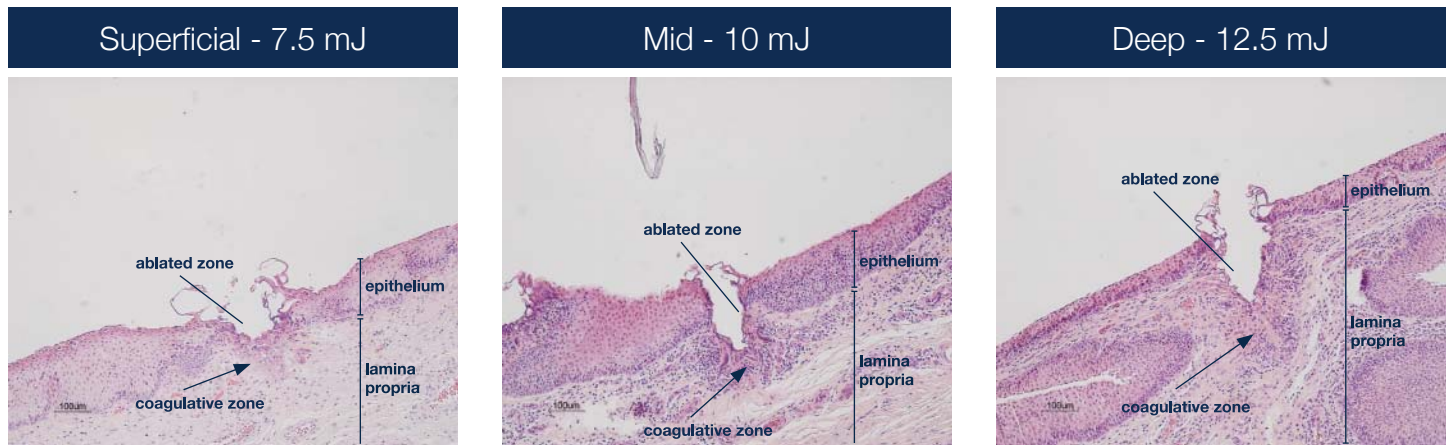

Histologies show total penetration depth (ablation and coagulation) at various energy levels of the AcuPulse treatments in the vaginal wall, using FemTouch.

* Ex-vivo porcine study

The AcuPulse product family is an upgradable comprehensive line of CO₂ Laser systems which your practice can grow with.

The wide range of accessories and clinical indications make the AcuPulse product line flexible and cost effective.

AcuPulse family of lasers extends your performance capabilities to treat a broader set of clinical applications, as well as more patients, in both the gynecology and aesthetic segments.

Gynecology

When it comes to operating on delicate anatomy, aiming for minimal thermal damage, the Lumenis CO₂ laser technology is perceived as superior amongst leading gynecologists around the world

Colposcopy and Vaginal approach

- › Cervical Intraepithelial Neoplasia (CIN)
- › Leukoplakia
- › Condyloma
- › Benign & Malignant tumors
- › Vaginal and vulvar resurfacing

Laparoscopy approach

- › Endometriosis
- › Adhesions
- › Myomas
- › Cysts
- › Hysterectomy

Aesthetics

AcuPulse has a complete ablative resurfacing workstation, featuring superior technology for all aesthetic resurfacing needs

- › Skin resurfacing
- › Acne scars
- › Dyschromia
- › Uneven texture
- › Actinic keratosis
- › lines and wrinkles
- › Uneven pigmentation
- › Skin furrows

Nearly **50%** of postmenopausal women suffer from stress urinary incontinence¹ and nearly **80%** suffer from vaginal atrophy²

"I am excited to add FemTouch treatment to my practice as this comfortable and easy procedure was demonstrating high satisfaction with my patients. This technology leads to good clinical outcomes and I can expand its utilization to additional gynecological indications"

Dr. Massimiliano Marziali, GYN, Rome, Italy



Let your patients

Laugh, Love, Live



FemTouch Specifications

Laser source	CO ₂
Wavelength	10,600 nm
Power mode	Continuous Wave (CW)
Laser delivery	Fractional
Pulse mode	Single, Repeat
Energy	7.5-12.5mJ
Density	3-25%
Certifications and approvals	Lumenis is ISO 13485:2012 certified, FemTouch on AcuPulse is CE approved and FDA cleared.
Compatible with	All AcuPulse and AcuPulse DUO models (except 40 WG)

¹ Melville JL. Urinary Incontinence in US Women, A Population-Based Study. Arch Intern Med. 2005. ² Simon JA. et al. Vaginal health in the United States: results from the vaginal health: insights, views & attitudes survey. Menopause (10) p1043-1048, 2013.

Manufactured by Lumenis Ltd.
Yokneam Industrial Park
6 Hakidma Street
P.O.B. 240
Yokneam 2069204, Israel
Tel: +972-4-959-9000



Lumenis (Germany) GmbH
Heinrich-Hertz-Str 3 D-63303
Dreieich-Dreieichenhain
GERMANY
Tel: +49 (0) 6103 8335 0

AMERICAS
San Jose, CA, USA
T +1 408 764 3000
+1 877 586 3647
F +1 408 764 3999

EMEA
Dreieich Dreieichenhain,
Germany
T +49 6103 8335 0
F +49 6103 8335 300

Roma (RM), Italy
T +39 06 90 75 230
F +39 06 90 75 269

Hertfordshire, UK
T +44 20 8736 4110
F +44 20 8736 4119

JAPAN
Tokyo, Japan
T +81 3 4431 8300
F +81 3 4431 8301

ASIA / PACIFIC
Beijing, China
T + 86 10 5737 6677

Gurgaon, India
T + 91 124 422 07 95

Kowloon, Hong Kong
T + 852 217 428 00
F + 852 272 251 51



WARNING: For safety use, Indications and Contraindications refer to the Operator Manuals and Instructions for Use of Lumenis Laser system and accessories.

© 2016 All Rights Reserved. The Lumenis Group of Companies. PB-2005136 rev A

www.Femtouch.com