

Pollogen  
BY ZUMENIS

DIVINE  
PRO

THE ULTIMATE FACIAL  
CONTOURING PLATFORM

# THE ULTIMATE MULTI-TREATMENT PLATFORM FOR FACIAL CONTOURING

The best RF in tri-powerful forms: Volumetric, Fractional, Resurfacing. One platform - all 3 layers of the skin.



## BENEFITS

### IMMEDIATE RESULTS THAT STAY

Effective treatments that provide immediate and long-lasting results, addressing the signs of aging by reducing wrinkles, tightening the skin and smoothing the face and the neck.

### COMFORTABLE AND ACCESSIBLE

Minimal pain, no downtime and doesn't require any anesthetics or mechanical force.

### SAFETY AND CLINICAL PROOF

Tested and validated technologies designed to stimulate the body's natural healing and reconstruction mechanisms.

### PERSONALIZED

A tailored treatment plan that addresses their needs and desires.

FOR ADDITIONAL CLINICAL DATA PLEASE VISIT [POLLOGEN.COM](https://www.pollogen.com)

# MARVEL AT THE ABILITIES OF OUR SKIN

Designed for comfort and flexibility, Divine Pro provides treatment combinations for encouraging and enhancing the body's own rejuvenation processes to achieve smoother, tighter skin, improved facial contouring and refreshed appearance of the skin.



## TRIPOLLAR RF

**Triggering regeneration of collagen and elastin fibers** by delivering volumetric RF energy safely and effectively, inducing controlled heating of the deep dermis.




RF Applicator



## VOLUDERM FRACTIONAL RF

Dual RF action, combined of a high energy pulse to penetrate the Stratum Corneum barrier, followed by a low energy pulse for ablation and coagulation of the viable underlying layers, resulting in **a focused, stronger effect in the dermis, with minimal damage to the epidermis and minimal pain or bruising.**

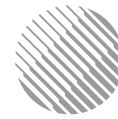

## TRIFRACTIONAL

Contact RF ablation for skin resurfacing and pin-point ablations of the epidermal layer, **inducing the natural wound healing process for skin remodeling and texture improvement**



Voluderm Applicator



GEN 12



GEN 36



GEN 36L



GEN 100



GEN 7X7

# DIVINE PRO ALTERNATING PROTOCOL

A highly effective combination of the Divine Pro technologies created for individuals with delicate, thin skin and/or sensitivities. Follow the two interchanging sessions, using the suggested technology combinations, with 3-4 weeks treatment-free intervals.

## Session A

<b>TECHNOLOGY</b>	<b>VOLUDERM</b>	<b>TRIPOLLAR RF</b>
<b>ENERGY TIP</b>	Dual Microneedling RF pulse via <b>Gen36</b> tip, with 6x6 array of Ultra-thin electrodes	Non-invasive Radiofrequency energy
<b>EFFECT</b>	Dermal coagulation	Skin tightening and volumizing

## Session B

<b>TECHNOLOGY</b>	<b>TRIFRACTIONAL</b>	<b>TRIPOLLAR RF</b>
<b>ENERGY TIP</b>	Single Resurfacing RF pulse via <b>H7X7</b> tip, with 7x7 array of electrode-pins	Non-invasive Radiofrequency energy
<b>EFFECT</b>	Intensive Epidermal ablation and coagulation with dermal coagulation	Skin tightening and volumizing

# DIVINE PRO LAYERING PROTOCOL

Follow this “Pyramid” treatment, layering all the available technologies, to treat all skin layers. For best results, it is recommended to move from the deepest skin layer towards the top layer: deep dermis, dermis, and epidermis. This treatment plan ensures the full exploitation of all the benefits offered by the Divine Pro platform.

Note: this plan is considered an advanced treatment for experienced users; therefore, it is advised to assess each patient by treating them with each technology separately before combining them all in a single session.

## Layering in same session

<b>TECHNOLOGY</b>	<b>TRIPOLLAR RF</b>	<b>VOLUDERM</b>	<b>TRIFRACTIONAL</b>
<b>ENERGY TIP</b>	Non-invasive Radiofrequency energy	Dual Microneedling RF pulse via <b>Gen100</b> tip, with 10x10 array of Ultra-thin electrodes or <b>Gen36L</b> tip, with 6x6 array of Ultra-thin electrodes	Single RF pulse via <b>H7X7</b> tip, with 7x7 array of electrode-pins
<b>EFFECT</b>	Skin tightening and volumizing	Dermal ablation and coagulation with additional epidermal effect	Intensive epidermal ablation and coagulation with additional dermal coagulation



## ABOUT POLLOGEN

Pollogen Ltd, a company of Lumenis BE, is a global leader in the medical and home-use aesthetics devices markets, providing innovative, safe, and effective solutions that enhance the practice of medical aesthetics professionals, as well as serving private aesthetic consumers. The company offers products that support a range of skin treatments, based on its proprietary technologies which include TriPollar, Volulderm and OxyGeneo.

Pollogen provides solutions for professionals and consumers in more than 60 countries, addressing the specific needs and the ever-changing standards of professionals and consumers in the dynamic beauty industry.

## ABOUT LUMENIS BE

Lumenis is a global leader in the field of minimally-invasive clinical solutions for the Surgical, Ophthalmology and Aesthetic markets, and a world renowned expert in developing and commercializing innovative energy-based technologies, including Laser, Intense Pulsed Light (IPL) and Radio-Frequency (RF).

### SYSTEM SPECIFICATIONS

**MAXIMUM OUTPUT RF POWER** Fractional Handpiece: ~50W TriPollar+® Handpiece: ~15W

**OUTPUT RF FREQUENCY** 1MHz

**INPUT VOLTAGE AND CURRENT** 100-240 Volt, 50-60 Hz, 1.4-0.7 A

**OUTPUT POWER RATING OF POWER SUPPLY** 24V DC, 5A

**USER INTERFACE** LCD TFT color touch screen

**DIMENSIONS** L32 x W32 x H43 cm; L12.6 x W12.6 x H17 in.

**WEIGHT** ~5 Kg; ~11 Lb


**HANDPIECES** Fractional Handpiece: RFF Microneedling & RFF Resurfacing  
TriPollar+® Handpiece: TriPollar®

**DISPOSABLE TIPS** Gen36 6x6 | Gen12 6x2 | H7x7 | Gen36L 6x6 | Gen100 10x10

© 2022 Pollogen Ltd. All rights reserved.

Trademarks and logos of Divine Pro, Tripollar, Volulderm and Pollogen are all trademarks of Pollogen Ltd.

Tripollar and Volulderm are proprietary technologies of Pollogen Ltd.

 Obelis s.a, Boulevard Général Wahis 53, 1030 Brussels, Belgium mail@obelis.net  
T +(32) 2.732.59.54 F +(32) 2.732.60.03

11600786 - rev.02 - March 2023